




# Supreme Gel

The Prism Supreme Gel™ Cushion was designed for users requiring mild to moderate positioning support and provides excellent pressure distribution relief. The gel overlay reduces shear and friction that contributes to skin breakdown, while the anatomically designed firm polyurethane base provides superior pelvic stability and leg positioning. A central core of pressure relieving visco-elastic insert under the IT's (Ischial Tuberosity), allows for excellent immersion for bony prominences.

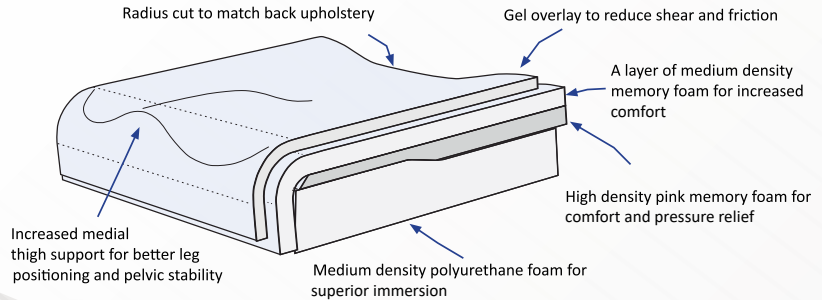

**Custom Sizes - Shipment within 7 business days**

**ADP CODE: SESFM0120**

Part Number	Size
SGCUS	Custom

**Non-Stock Sizes – Shipment within 3 business days**

**ADP CODE: SESFM0110**

Part Number	Size	Part Number	Size
SGC1516	15" x 16"	SGC1916	19" x 16"
SGC1518	15" x 18"	SGC1918	19" x 18"
SGC1620	16" x 20"	SGC1920	19" x 20"
SGC1716	17" x 16"	SGC2016	20" x 16"
SGC1718	17" x 18"	SGC2018	20" x 18"
SGC1720	17" x 20"	SGC2020	20" x 20"

**Stock Sizes – Same Day Shipment**

**ADP CODE: SESFM0115**

Part Numb	Size	Part Number	Size
SGC1616	16" x 16"	SGC1818	18" x 18"
SGC1618	16" x 18"	SGC1820	18" x 20"
SGC1816	18" x 16"		

## Product Highlights

- Gel overlay provides even distribution of weight across the cushion
- Waterfall front edge for foot propelling comfort
- Cane cutout in rear allowing correct placement in wheelchair and proper seat depth
- 3" of visco foam under IT's for superior immersion
- Supreme pressure relief
- Non slip cover
- Standard 4-way stretch, water resistant inner liner
- Standard 4-way stretch, breathable outer cover

Product Weight: 10 lbs



User Weight: 250 lbs

Warranty: Two (2) years limited on cushion and 90 days on the cover.


### Best Used For:

Adv. Weight Distribution  
Heat Management  
Extended Length Use  
Adv. Positioning Needs  
Maintain Skin Integrity  
Reduce Fluid Migration


Pressure Reading Scale in mmHG

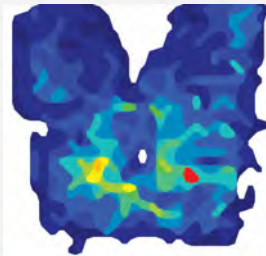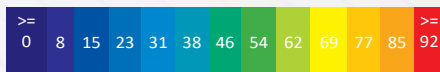

Pressure Map	Supreme Gel
Weight	150 Lbs/68 Kg
Duration	1 hr
Average Pressure	24 mmHg

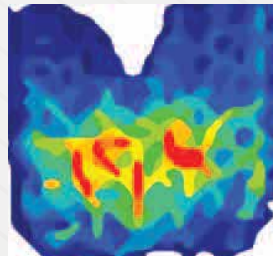


Pressure Map	Supreme Gel
Weight	220 Lbs/99 Kg
Duration	1 hr
Average Pressure	35 mmHg

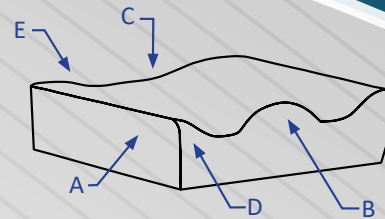

### CUSHION THICKNESS

- A- Height at adductor = 5"
- B- Height at abductor = 5"
- C- Height at rear IT = 3<sup>1/4</sup>"
- D- Leg Trough = 3<sup>1/4</sup>"
- E- Height at rear trochantar = 4<sup>1/4</sup>"

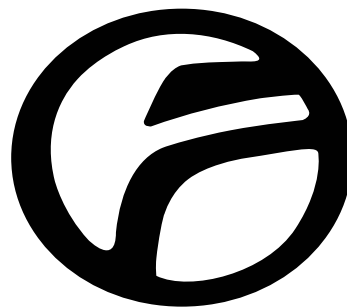

**FUTURE MOBILITY**

**CANADA**

6750 Professional Court, Mississauga, L4V1X6

1-888-737-4011

[www.futuremobility.ca](http://www.futuremobility.ca)

**USA**

1 Buffalo River Place, Buffalo, NY 14210

1-855-938-8873

[www.futuremobility.com](http://www.futuremobility.com)