




Supreme II

The Prism Supreme™ was designed for users requiring moderate to high positioning support and provides excellent pressure distribution relief. The anatomically designed firm polyurethane base, with its aggressive leg troughs and abductor contouring, provides superior pelvic stability and leg positioning. A central core of pressure relieving visco-elastic insert under the IT's (Ischial Tuberosity), allows for excellent immersion for bony prominences.

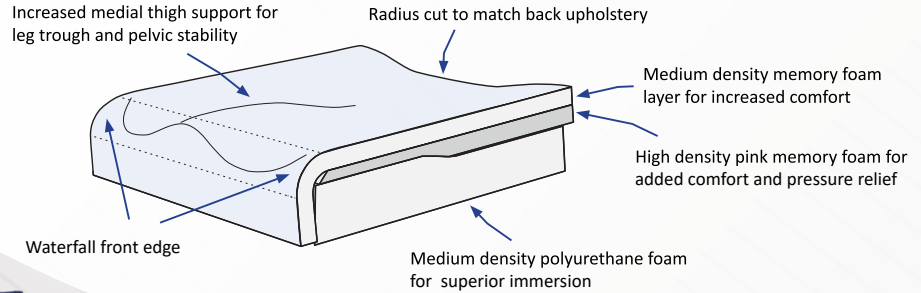

Custom Sizes - Shipment within 7 business days

ADP CODE: SESFM0025

Part Number	Size
SCCUS	Custom

Non-Stock Sizes – Shipment within 3 business days

ADP CODE: SESFM0020

Part Number	Size	Part Number	Size
SC1516	15" x 16"	SC1916	19" x 16"
SC1518	15" x 18"	SC1918	19" x 18"
SC1620	16" x 20"	SC1920	19" x 20"
SC1716	17" x 16"	SC2016	20" x 16"
SC1718	17" x 18"	SC2018	20" x 18"
SC1720	17" x 20"	SC2020	20" x 20"

Stock Sizes – Same Day Shipment

ADP CODE: SESFM0015

Part Number	Size	Part Number	Size
SC1414	14" x 14"	SC1618	16" x 18"
SC1416	14" x 16"	SC1818	18" x 18"
SC1616	16" x 16"	SC1820	18" x 20"

Product Highlights

- 3" of visco foam (memory foam) under IT's (Ischial Tuberosity) for superior immersion
- High density, contour cut layered foam for a consistent, long lasting product
- 4-way stretch, water resistant inner liner
- Radius cut to match shape of back upholstery
- 4-way stretch, breathable outer cover

- Waterfall front edge for foot propelling comfort
- Cane cutout in rear to allow for correct placement in wheelchair and proper seat dept

Product Weight: 3.5 lbs

User Weight: 250 lbs

Warranty: Two (2) years limited on cushion and 90 days on the cover.

Supreme Modifications

Part Number	Size	Part Number	Size
SLDL1	Left 1"	SLDR1	Right 1"
SLDL2	Left 2"	SLDR2	Right 2"



Pelvic Obliquity Build-Up

Part Number	Size	Part Number	Size
POBL.5	Left 1/2"	POBR.5	Right 1/2"
POBL1	Left 1"	POBR1	Right 1"
POBL1.5	Left-1-1/2"	POBR1.5	Right 1-1/2"


Best Used For:

Adv. Weight Distribution
Heat Management
Extended Length Use
Adv. Positioning Needs
Maintain Skin Integrity
Asymmetrical Posture


Pressure Reading Scale in mmHG

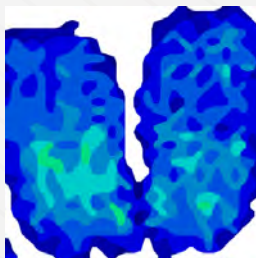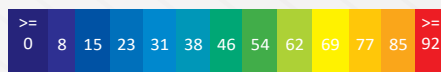

Pressure Map	Supreme II
Weight	125 Lbs/57 Kg
Duration	1 hr
Average Pressure	25 mmHg

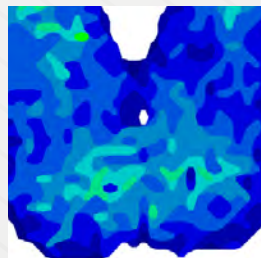


Pressure Map	Supreme II
Weight	200 Lbs/91 Kg
Duration	1 hr
Average Pressure	30 mmHg

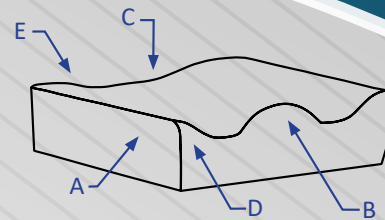

CUSHION THICKNESS

- A- Height at adductor = 5"
- B- Height at abductor = 5"
- C- Height at rear IT = 3^{1/4}"
- D- Leg Trough = 3^{1/4}"
- E- Height at rear trochanter = 4^{1/4}"

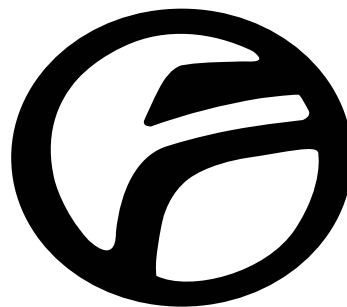

FUTURE MOBILITY

CANADA

6750 Professional Court, Mississauga, L4V1X6

1-888-737-4011

www.futuremobility.ca

USA

1 Buffalo River Place, Buffalo, NY 14210

1-855-938-8873

www.futuremobility.com