


Super Cushion For Kids

The **Super Cushion for Kids™** is a Junior size cushion created with three layers of polyurethane, high density and memory foam, which provides children with excellent pressure relief. The anatomically designed firm foam base with its aggressive leg trough and abductor contour, provides superior pelvic stability and leg positioning. A central core of pressure relieving visco-elastic insert under the IT's (Ischial Tuberosity), allows for excellent immersion for bony prominences.

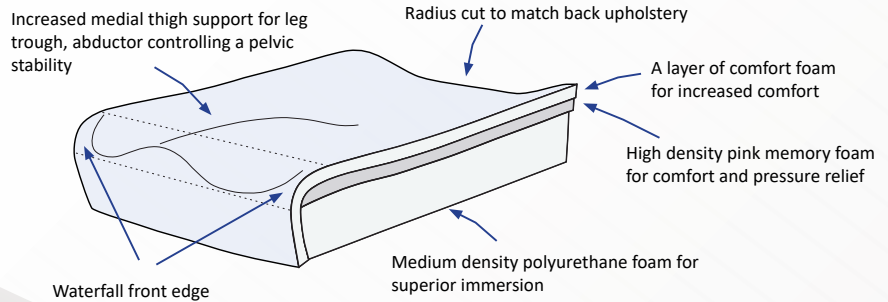

Custom Sizes - Shipment within 7 business days

ADP CODE: SESFM0100

Part Number	Size
KSCCUS	Custom

Non-Stock Sizes – Shipment within 3 business days

ADP CODE: SESFM0095

Part Number	Size	Part Number	Size
KSC1012	10" x 12"	KSC1216	12" x 16"
KSC1014	10" x 14"	KSC1410	14" x 10"
KSC1016	10" x 16"	KSC1412	14" x 12"
KSC1210	12" x 10"	KSC1610	16" x 10"
KSC1214	12" x 14"	KSC1612	16" x 12"

Stock Sizes – Same Day Shipment

ADP CODE: SESFM0090

Part Number	Size	Part Number	Size
KSC1010	10" x 10"	KSC1416	14" x 16"
KSC1212	12" x 12"	KSC1614	16" x 14"
KSC1414	14" x 14"	KSC1616	16" x 16"

Product Highlights

- High density layered foam for a consistent, long lasting product
- Standard with water resistant outer cover, made of black Meshtex urethane
- Waterfall front edge for foot propelling comfort
- Radius cut to match shape of back upholstery
- Pressure relief
- Optional 4-way stretch, water resistant inner liner
- 4-way stretch, breathable outer cover


Product Weight: 1.5 lbs

User Max Weight: 150 lbs

Warranty: Two (2) years limited on cushion and 90 days on the cover.

Best Used For:



Pediatric Specialized
Cerebral Palsy
Postural Asymmetries
Adv. Positioning Needs
Muscular Dystrophy
Maintain Skin Integrity

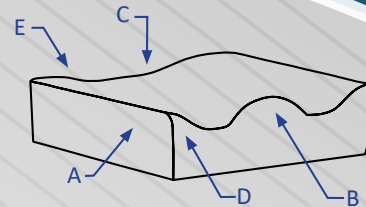


CUSHION THICKNESS

- A- Height at adductor = 3 3/4"
- B- Height at abductor = 4"
- C- Height at rear IT = 2 1/2"
- D- Leg Trough = 3 1/2"
- E- Height at rear trochanter = 3 1/2"

Pressure Reading Scale in mmHG

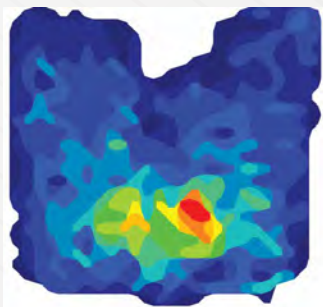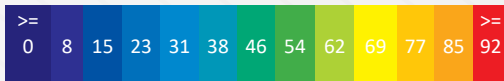

Pressure Map	Super Kids
Weight	150 Lbs/68 Kg
Duration	1 hr
Average Pressure	24 mmHg

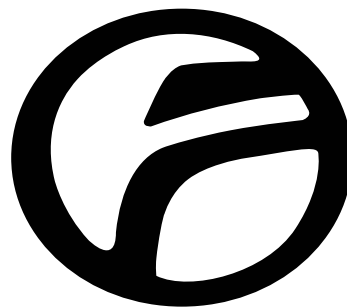

FUTURE MOBILITY

CANADA

6750 Professional Court, Mississauga, L4V1X6

1-888-737-4011

www.futuremobility.ca

USA

11126 Greenstone Ave, Santa Fe Springs, CA 90670

716-783-9130

www.futuremobility.com